


Chromoplasts

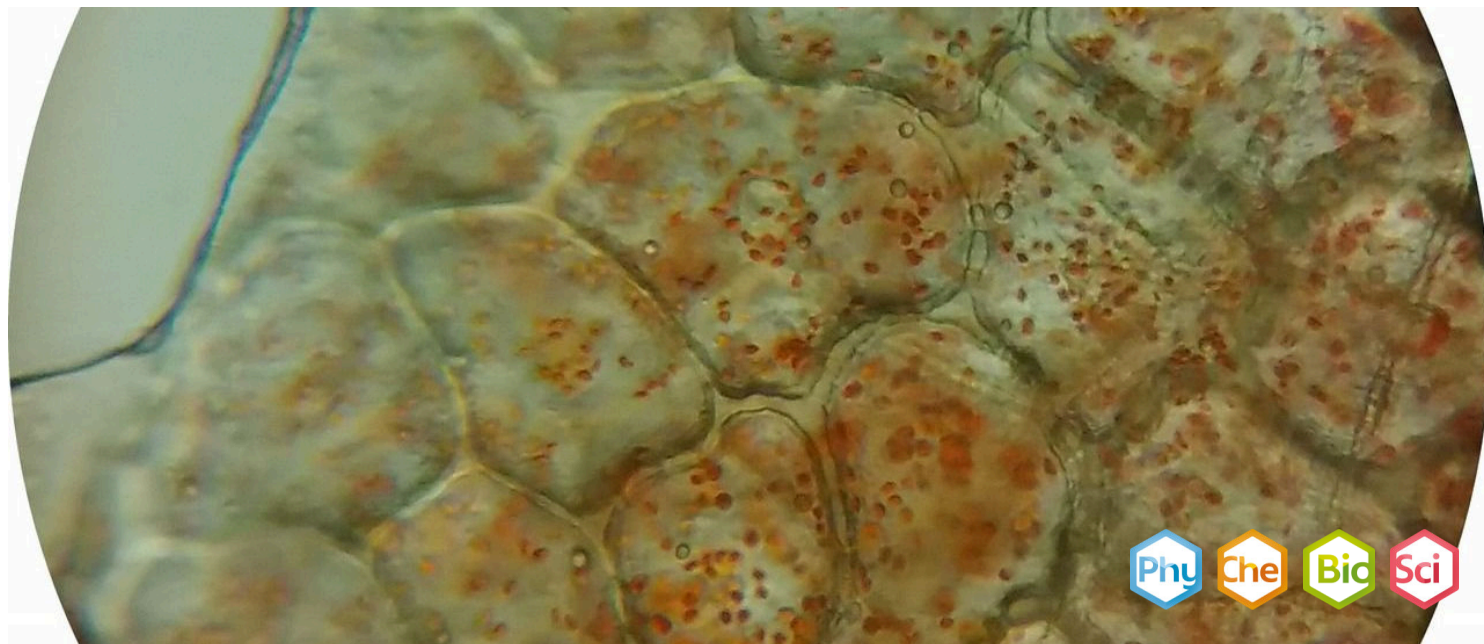


Biology → Microscopy / Cell Biology → Basics of Microscopy & Work Technology

Biology → Microscopy / Cell Biology → Plants & Fungi

Biology → Microscopy / Cell Biology → Cell structure

Nature & technology → From the very small & the very big

Nature & technology → Plants & animals


Difficulty level

easy


Group size

1


Preparation time

10 minutes


Execution time

30 minutes

This content can also be found online at:



<http://localhost:1337/c/5f5126e5739d0a0003ee3e74>


PHYWE



Teacher information

Application

PHYWE

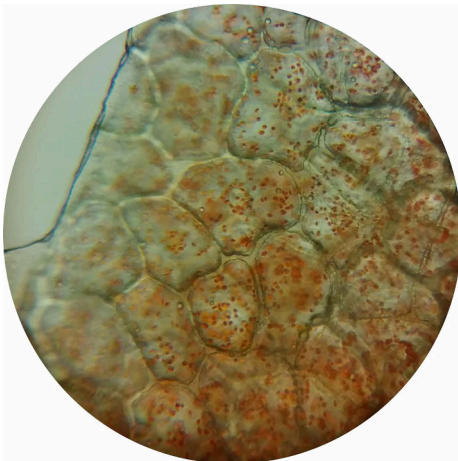

Cells of peppers at low magnification

The fruits and blossoms of many plants have strong, bright colours. The coloured fruits attract animals, are eaten by them and the seeds are excreted elsewhere. This is how the plant is spread. Coloured flowers attract insects, which harvest the nectar from the flowers. In addition, the pollen is transported from plant to plant, thus ensuring fertilization. Yellow and red dyes are usually found in certain cell organelles, the chromoplasts.

Other teacher information (1/3)

PHYWE

Prior knowledge



The cytoplasm of plant cells contains numerous plastids surrounded by thin membranes. They are usually ovoid in shape and are distinguished by their colour. Due to photosynthesis, the most important group is that of chloroplasts. The chromoplasts contain carotenes or xanthophylls. Some of the colours have a meaning for attracting animals (e.g. flowers - insects, fruits - birds).

Scientific principle



The students should understand the preparation instructions as a stimulus. They can make peel preparations, cut preparations, but also scraping preparations of carrots or squeezing preparations of the soft flesh of tomatoes. It is important to make sure that the preparations become thin, i.e. that they are translucent.

Other teacher information (2/3)

PHYWE

Learning objective



The students learn how chromoplasts are present in form and arrangement in the cell.



Tasks



Making and microscoping the preparation

1. Red pepper - fruit skin
2. Capuchin Cress - Petal
3. Rose - pulp

Other teacher information (3/3)

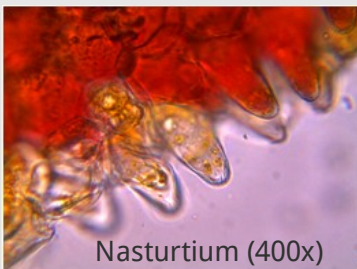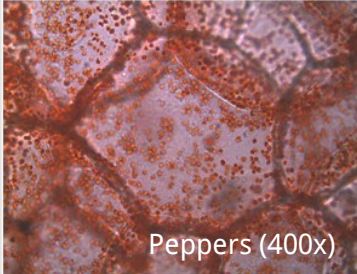

Chromoplasts are present in all yellow, red and orange parts of the plant. The preparations suggested on the student worksheet are particularly suitable, but students should be encouraged to examine fruits and flowers that are easily accessible to them. The pupils can bring coloured vegetables from the kitchen at home and flowers from the garden. Blue and purple flowers are not suitable, as the most common dye anthocyanin is dissolved in the vacuoles. Fruits of the rose can be harvested from wild roses or from bed roses.

Safety instructions

PHYWE



- Razor blades are sharp. Extreme caution is advised when handling them.
- Working with microscopes for too long can lead to physical discomfort (fatigue, headaches, nausea), especially when students are inexperienced.
- Microscopes are sensitive. During transport and handling, care should be taken to ensure that everything is done carefully and without rushing.
- The general instructions for safe experimentation in science teaching apply to this experiment.


PHYWE



Student Information

Motivation

PHYWE

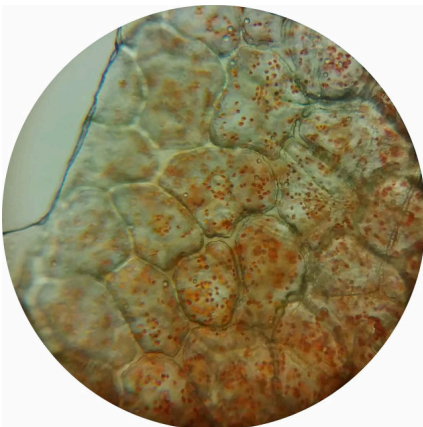


Cells of peppers at low magnification


In this experiment you will learn how to recognize the chromoplasts in different preparations of different plants and fruits and how to describe their shape and arrangement.

Tasks

PHYWE

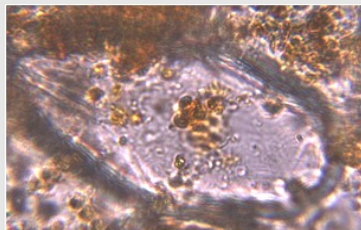

Fruits of the rose (400x)

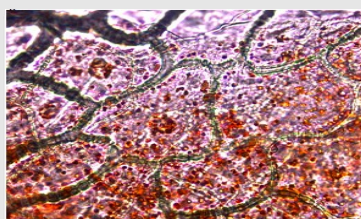

Peppers (400x)

Preparation and microscopy of preparations


1. Red pepper - fruit skin
2. Capuchin Cress - Petal
3. Rose - pulp

Equipment

Position	Material	Item No.	Quantity
1	PHYWE Binocular student microscope, 1000x, mechanical stage	MIC-129A	1
2	Microscopic slides, 50 pcs	64691-00	1
3	Cover glasses 18x18 mm, 50 pcs	64685-00	1
4	Beaker, 100 ml, plastic (PP)	36011-01	1
5	Dropping pipette with bulb, 10pcs	47131-01	1
6	Tweezers, straight, pointed, 120mm	64607-00	1
7	Scalpel holder	64615-00	1
8	Scalpel blades, rounded tip, 10 off	64615-02	1

Procedure (1/3)

PHYWE

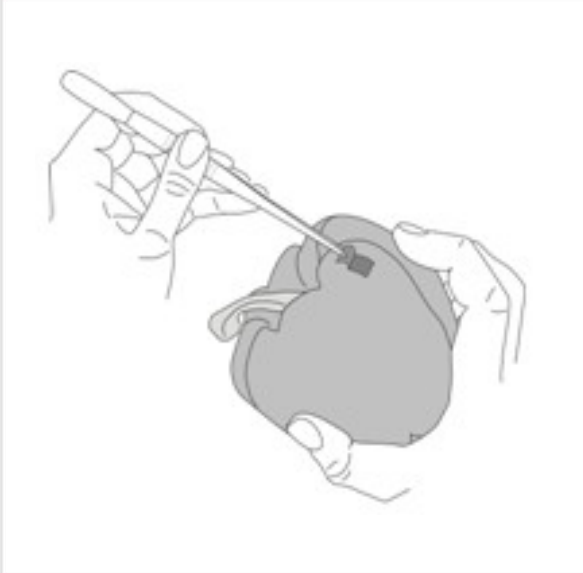

(1) Prepare the sample: Red pepper - fruit skin

- Prepare a slide with a drop of water.
- With the tweezers a piece of the skin of the pepper fruit is peeled off.
- The cuticle is placed directly into the water drop, covered and microscopically examined.

Procedure (2/3)


(2) Prepare the sample: Capuchin cress - petal

- Prepare a slide with a drop of water.
- A yellow or red petal is placed over the finger and a thin surface incision is made with the scalpel.
- The preparation is placed directly into the water drop, covered and microscopically examined.



Procedure (3/3)

PHYWE

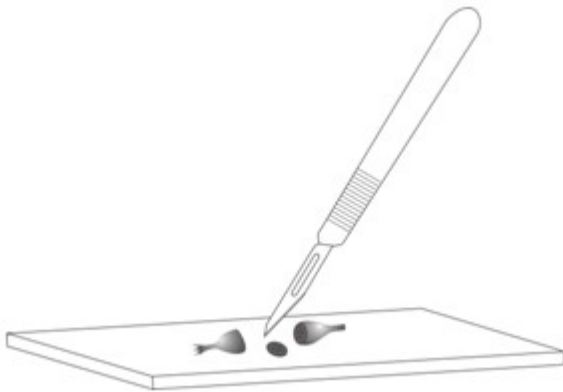


(3) Prepare the sample: Rose - pulp

- Prepare a slide with a drop of water.
- From the red fruit of a rose (rosehip) a very thin piece of fruit is cut off. The even flesh of the fruit is examined, not the seed.
- The preparation is placed directly into the water drop, covered with a cover glass and microscopically examined.

PHYWE

Report

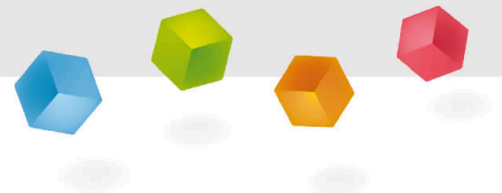

Task 1


PHYWE

In the plant parts with which colours are chromoplasts contained?

☐ Greens☐ Yellow☐ Oranges☐ Red☒ Check

Task 2

PHYWE

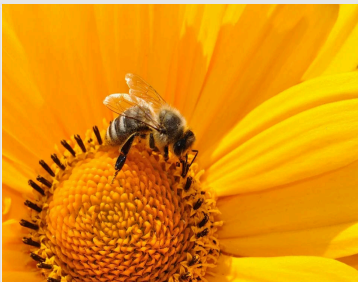


What is a possible meaning of the colours in the chromoplasts?

Task 3

PHYWE


Microscopy with increasing magnification. Drawing from the specimen **Red pepper - fruit skin** a cell. Use a coloured pencil to show the location and number of chromoplasts.

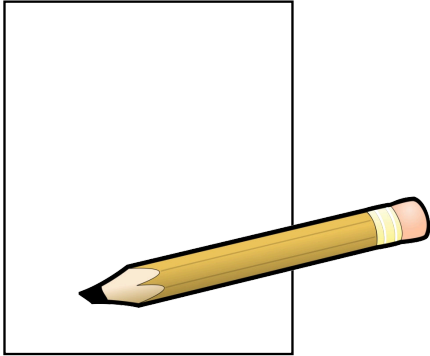

Task 4

PHYWE


Microscopy with increasing magnification. Drawing from the specimen **Capuchin Cress - Petal** a cell. Use a coloured pencil to show the location and number of chromoplasts.

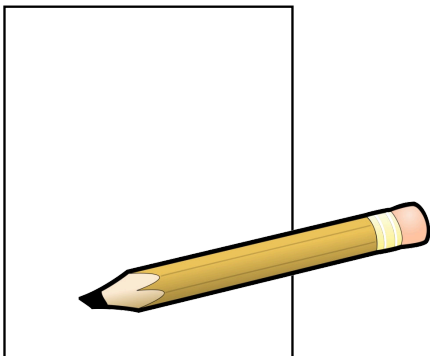

Task 5

PHYWE

Microscopy with increasing magnification. Drawing from the specimen **Rose - pulp** a cell. Use a coloured pencil to show the location and number of chromoplasts.

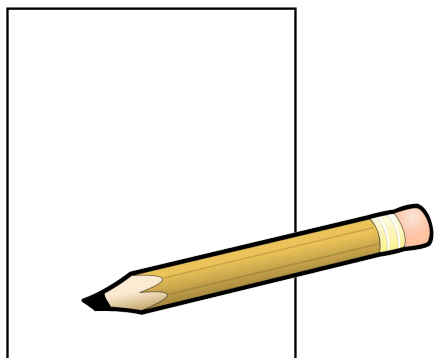

Slide

Score/Total

Slide 15: Plant parts with chromoplasts

0/3

Slide 16: Colours in chromoplastics

0/1

Total amount

 0/4 Solutions Repeat