

Cooling by evacuation with SMARTsense



Students will be able to understand the cooling by evacuation process.

Physics

Thermodynamics

Temperature & Heat



Difficulty level

medium



Group size

2



Preparation time

10 minutes



Execution time

20 minutes


This content can also be found online at:



<http://localhost:1337/c/60b220dd6b5a0f00039aff31>

PHYWE

General information

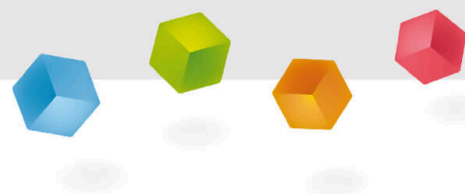

Application

PHYWE



Experimental setup

In the laboratory, water evaporates at an air pressure of 1.013 bar up to a temperature of 100 °C; only above 100 °C water boils. Evaporation is determined by the temperature of the liquid and air pressure.

Evaporation requires energy, which is mainly taken from the evaporated substance. The evaporating liquid cools down in the process. In general: the lower the air pressure, the lower the boiling temperature and the higher the evaporation rate.

In this experiment students examine the correlation between air pressure and temperature of a liquid. They determine the temperature curve of water during lowering the air pressure by a vacuum pump.

Other information (1/2)

PHYWE

Prior knowledge



Students should already be familiar with the principles of vacuum, gas pressure and boiling temperatures.

Scientific principle



When the air pressure above a water surface is reduced, the water begins to boil at a certain temperature. The temperature of the water is hereby reduced and further evacuation can finally bring it to 0 °C and even lower.

Other information (2/2)

PHYWE

Learning objective



Students will be able to understand the cooling by evacuation process. They learn about the influence of the air pressure on a liquid (and their temperature). Students can explain the cooling effect with the kinetic gas theory.

Tasks



The students will determine the temperature curve of water during pumping.

Safety instructions

PHYWE



- For this experiment the general instructions for safe experimentation in science lessons apply.

Theory

PHYWE


When a liquid boils, molecules of the liquid exit from the liquid into the gas space above it. The heat of evaporation that is required for this is drawn from the water. To be able to exit from the water, the molecules of the liquid must have a defined kinetic energy that is greater the higher the gas pressure above the liquid is.

According to the kinetic gas theory, the square of the average speed of the molecules of the liquid is proportional to the temperature of the liquid. This means on the one hand that the boiling temperature is reduced at lower gas pressures, and on the other hand that the molecules having the most kinetic energy exit the liquid, so that the temperature of the remaining liquid is reduced. The "hottest" water molecules are pumped off, so that average speed of the molecules of liquid, and so the temperature of the liquid, continually decreases.

Equipment

Position	Material	Item No.	Quantity
1	Cobra SMARTsense - Temperature, - 40 ... 120 °C (Bluetooth)	12903-00	1
2	Rubber tubing,vacuum,i.d.8mm	39288-00	2
3	Hose clamp for 12-22 mm diameter	41002-00	2
4	Rotary valve vacuum pump, two stages, 115 V / 230 V	02741-95	1
5	Pump plate, complete	02668-88	1
6	Support base DEMO	02007-55	1
7	Bell jar, with knob and sealing ring	02668-10	1
8	Protection cylinder for bell-jar	02668-14	1
9	Beaker, Borosilicate, tall form, 1000 ml	46030-00	1
10	Boiling beads, 200 g	36937-20	1
11	measureAPP - the free measurement software for all devices and operating systems	14581-61	1
12	Felt sheet, 100 x 100 mm	04404-20	2

PHYWE

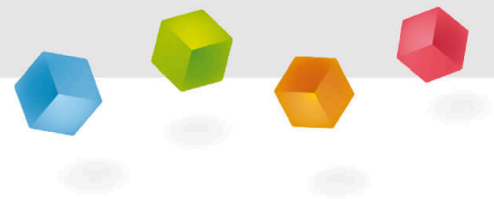

Setup and procedure

Setup (1/3)

PHYWE

To measure the temperature, you need the Cobra SMARTsense Temperature sensor and the measureAPP. Check whether "Bluetooth" is activated on your device (tablet, smartphone) (the app can be downloaded free of charge from the App Store - QR codes below). Now open the measureAPP on your device.



measureAPP for

Android operating systems



measureAPP for

iOS operating systems




measureAPP and

Tablets and PCs with Windows 10

Setup (2/3)

PHYWE

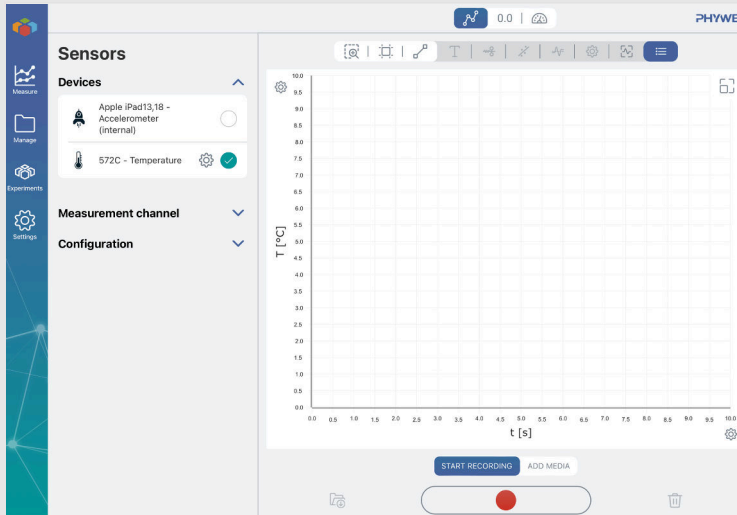

Sensor "Temperature"

- Switch on the SMARTsense sensor "Temperature" by pressing and holding the power button.
- Connect the sensor to the device in the measureAPP under the item "Measure" as shown in the illustration on the left.
- The SMARTSense sensor "Temperature" is now displayed in the app.

Setup (3/3)

PHYWE




Experimental setup

- Set up the experiment as shown in Fig. right.
- Fill about 100 ml of water (room temperature) in the beaker and add two boiling stones.
- Place the beaker on the pump plate and lay the felt sheets under the beaker as insulation.
- Dip the temperature sensor in the water.
- Place the bell jar with rubber sealing ring and the protection cylinder on the pump plate.
- Fit the oil mist filter to the vacuum pump and connect the vacuum pump to the manometer and the pump plate.

Procedure


PHYWE

- Start measuring the temperature with clicking on "Start".
- Switch on the vacuum pump and check on the manometer that there is no leak in the vacuum line that would be detrimental to evacuation. The pressure should drop below 100 mbar within the first minute of pumping.
- Switch the vacuum pump off after about 40 minutes.
- End measurement in the app with clicking on "Stop".

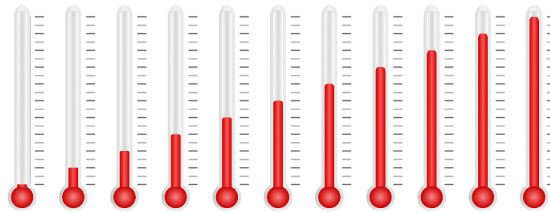

PHYWE

Evaluation

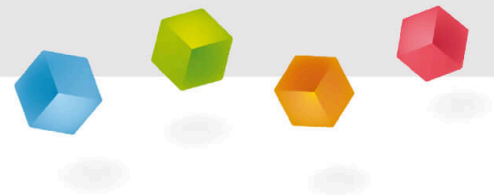

Evaluation (1/6)

PHYWE

Results

- It can be seen that the water starts to boil after about 1-2 minutes.
- After about five minutes the water has cooled by about 10°C.
- After about half an hour the temperature of the water reaches the freezing point, a layer of ice starts to form on the water surface.

When a liquid boils, molecules of the liquid exit from the liquid into the gas space above it. The heat of evaporation that is required for this is drawn from the water. To be able to exit from the water, the molecules of the liquid must have a defined kinetic energy that is greater the higher the gas pressure above the liquid is.

Evaluation (2/6)

PHYWE


According to the kinetic gas theory, the square of the average speed of the molecules of the liquid is proportional to the temperature of the liquid.

This means on the one hand that the boiling temperature is reduced at lower gas pressures, and on the other hand that the molecules having the most kinetic energy exit the liquid, so that the temperature of the remaining liquid is reduced. The "hottest" water molecules are pumped off, so that average speed of the molecules of liquid, and so the temperature of the liquid, continually decreases.

In the last part of the curve (Fig. 2, next slide), a slight rise in the measured temperature can be seen after about 1,800 seconds. The explanation for this is that the water is first supercooled a little, but then the surface froze whereby heat of fusion was released that led to the observed temperature rise measured at the bottom of the vessel.

Evaluation (3/6)

PHYWE

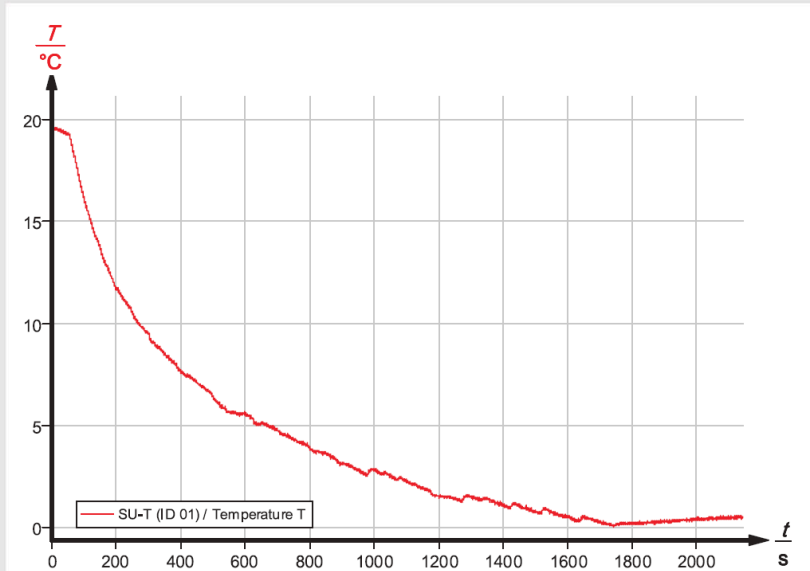

Temperature curve of the water during pumping. The boiling temperature of water is dependent on the pressure of the air that burdens it.

When the air pressure is reduced, the boiling point is also reduced. The boiling water then cools until it finally reaches a point on the steam pressure curve where it is again in equilibrium with the prevailing pressure. Cooling can so finally bring the temperature to below freezing point.

Evaluation (4/6)

PHYWE

What happens about half an hour after the experiment is started?

- ☐ After about half an hour the temperature of the water reaches the boiling point, a layer of steam starts to form on the water surface.
- ☐ After about half an hour the temperature of the water reaches room temperature.
- ☐ After about half an hour the temperature of the water reaches the freezing point, a layer of ice starts to form on the water surface.

✓ Überprüfen

Evaluation (5/6)

PHYWE

What exactly happens, when a liquid boils?

- ☐ When a liquid boils, molecules of the liquid exit from the liquid into the gas space above it.
- ☐ When a liquid boils, molecules of the liquid exit from the liquid into the ice phase on the surface.
- ☐ When a liquid boils, molecules of the liquid will stop moving around like they do when it freezes.
- ☐ When a liquid boils, nothing happens except a faster movement of the molecules.

✓ Überprüfen

Evaluation (6/6)

PHYWE


Choose the correct answers.

- ☐ The heat of evaporation that is required for the liquid to exit from liquid to gas space is drawn from the water.
- ☐ According to the kinetic gas theory, the square of the average speed of the molecules of the liquid is un-proportional to the temperature of the liquid.
- ☐ According to the kinetic gas theory, the square of the average speed of the molecules of the liquid is proportional to the temperature of the liquid.

✓ Überprüfen

Slide	Score / Total
Slide 17: After half an hour	0/1
Slide 18: Liquid boiling	0/1
Slide 19: Temperature of the liquid	0/2

Total Score

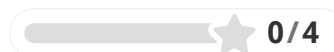

0/4



Show solutions



Retry