

PHYWE Systeme GmbH & Co. KG
Robert-Bosch-Breite 10
D-37079 Göttingen

Telefon +49 (0) 551 604-0
Fax +49 (0) 551 604-107
E-mail info@phywe.de
Internet www.phywe.com

Operating instructions


 The unit complies with the applicable EC-guidelines



Fig. 1: 12902-02 Cobra SMARTsense Current ± 1 A

TABLE OF CONTENTS

1	SAFETY PRECAUTIONS
2	PURPOSE AND CHARACTERISTICS
3	FUNCTIONAL AND OPERATING ELEMENTS
4	NOTES ON OPERATION
5	HANDLING
6	TECHNICAL DATA
7	SCOPE OF DELIVERY
8	ACCESSORIES
9	CONFORMITY
10	DISPOSAL
11	NOTES ON BATTERY AND RECHARGEABLE BATTERY DISPOSAL

1 SAFETY PRECAUTIONS



Caution!

- Carefully read these operating instructions completely before operating this instrument. This is necessary to avoid damage to it, as well as for user-safety.
- Only use the instrument for the purpose for which it was designed.
- Only use the instrument in dry rooms in which there is no risk of explosion.
- Protect the instrument from dust, moisture and vapours. Use a slightly moist lint-free cloth to clean the instrument. Do not use aggressive cleaning agents or solvents.
- Take care that no liquid penetrates in through the housing openings, as such penetration would result in damage to Sensor.
- Do not open the unit.
- Do not use the device in electrical circuits where voltage > 25 VAC or > 60 VDC may occur.
- The sensor is protected by way of a self-resetting fuse (a semiconductor element which, in the event of a high current load, becomes highly resistive, thereby reducing the current flow. The electrical circuit, however, is not completely interrupted!). Avoid short-circuit current above 20 A. It may destroy the sensor.

2 PURPOSE AND CHARACTERISTICS

The sensor is used for measuring current values and for transferring the values to a terminal device, e.g. a tablet computer, smartphone, etc., via Bluetooth or USB.

3 FUNCTIONAL AND OPERATING ELEMENTS


3.1 Operating elements

The sensor has an on-button and two LEDs for indicating the Bluetooth and battery charge status.

On-Button 

Pressed for longer 3s	Switch sensor on/off
Pressed 3x quickly	Start offline measurement
Pressed 2x quickly	Stop offline measurement

If the sensor is to be connected via USB, it is not necessary to press the power button longer 3s.

Bluetooth-LED 

Flashing red every 2 seconds	Not connected
Flashing green every 2 seconds	Connected to the terminal device
Flashing green every 4 seconds	Running measurement

Battery charge LED 

Flashing red every 2 seconds	Low battery
Illuminated red	Active charging process
Illuminated green	Charging process completed

3.2 Measurement inputs

The safety sockets for the current measurement are located on the right and left on the face of the sensor.



Due to the low internal resistance of the sensor, a short circuit may quickly lead to high current flowing through the sensor. If this current exceeds 20 A, even briefly, the sensor may be destroyed. This is why a correct experiment set-up is of utmost importance.

3.3 USB port

The battery, which is permanently installed in the sensor, is charged via the type C USB port. Furthermore, communication with the terminal device (Computer/Tablet) takes place via this interface.

4 NOTES ON OPERATION

The device fulfils all of the technical requirements that are compiled in current EC guidelines. The characteristics of this product qualify it for the CE mark.

This instrument is only to be put into operation under specialist supervision in a controlled electromagnetic environment in research, educational and training facilities (schools, universities, institutes and laboratories).

The individual connecting leads are each not to be longer than 2 m.

The instrument can be so influenced by electrostatic charges and other electromagnetic phenomena (HF, bursts, indirect lightning discharges) that it no longer works within the given specifications. Carry out the following measures to reduce or eliminate the effect of such disturbance: Ensure potential equalization at the PC (especially with Laptops). Use screening.

5 HANDLING

This section describes the start-up of the sensor and the recording of measurement data. Please read this section thoroughly in order to avoid failures or operating errors.

5.1 Charging process

Use a USB-C cable to connect the sensor to a computer or USB charger (not included).

During the charging process, the battery charge LED lights up red. When the charging process is complete, the battery charge LED lights up green. The charging time for a completely discharged battery is 3 hours maximum.



Disconnect the charger at the latest four hours after the completion of the charging process. Otherwise, the service life of the battery may be negatively affected.

5.2 Start-up

Switch on the sensor by pressing the power button for more than 3s. Now the Bluetooth LED flashes red. Start the software and select the sensor.

If the sensor is to be used via the USB interface, it does not need to be switched on. The sensor is connected directly to the end device using the supplied USB cable.

A 9-digit code is printed on the back of the sensor (Fig.2). The last 4 digits of the code are displayed in the software as the sensor designation (Fig.3). This allows an exact assignment of the sensors possible with the software.



Fig. 2

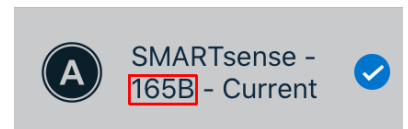


Fig. 3

Selection of the sensor via the Bluetooth interface

Make sure that the Bluetooth interface is activated on the terminal device (PC/Tablet/Smartphone) and that the software is allowed to access the interface.

After the sensor has been selected in the software, the LED flashes green to indicate that the connection has been established correctly. After the sensor has been coupled with the software, the sensor is no longer visible to other users in the software, and therefore can no longer be selected.

If the sensor is switched on and not connected, it switches off automatically after 5 minutes.

Selection of the sensor via the USB interface

For this purpose, the sensor must be plugged into the USB port of the end device. It is not necessary to switch on the sensor. The sensor is automatically recognized and displayed. It can be selected and connected directly.

5.3 Offline measurement

To start an offline measurement, press the power button 3 times in quick succession. The Bluetooth LED then flashes green 3 times in rapid succession to acknowledge the successful start. To stop a measurement, press the switch-on button 2x in quick succession. The Bluetooth LED also acknowledges this by flashing quickly.

Offline measurements can be read out via the measureAPP or measureLAB software. Furthermore, offline parameters such as data rate and measurement duration can be set. After the set measurement duration has elapsed, the offline measurement is automatically terminated. However, the measurement can always be ended prematurely by pressing the switch-on button.

6 TECHNICAL DATA

Operating temperature range: 5 - 40°C

Rel. humidity < 80%

Measuring range	±1 A
Resolution	0.5 mA
Accuracy*	± 1 %
Ri (typical)	300 mΩ
Max. data rate	100 kHz
Max. perm. Voltage in the circuit	25 V
Overcurrent protection	1.1A, self-resetting
Battery	LiPo 3,7V / 250 mAh
Max. wireless range (open field))	30 m
Dimensions (length x width x height)	90 x 55 x 23 mm
Weight	49 g

*referred to the full-scale value

7 SCOPE OF DELIVERY

- The extent of delivery is as follows
- Cobra SMARTsense Current 12902-02
- USB connecting cable type C 07935-00
- Operating instructions

8 ACCESSORIES

The following accessories are available:

- Cobra SMARTlink 12999-99
- Connecting cable, 32 A, 750 mm, red 07362-01
- Connecting cable, 32 A, 750 mm, blue 07362-04
- USB-charger 07934-99
- USB connecting cable type C 07935-00
- Software measureLAB 14580-61
- Free measureApp available from supplier portals

iOS



Android



Windows



9 CONFORMITY



PHYWE Systeme GmbH & Co.KG herewith declares that the radio system type 12902-02 complies with Directive 2014/53/EU. The complete text of the EU Declaration of Conformity is available at the following Internet address:

www.phywe.com/en/ec-declaration

10 DISPOSAL

The packaging consists mainly of environmentally friendly materials that should be disposed of at local recycling centres.



Do not dispose of this product with normal household waste. If this unit needs to be disposed of, please return it to the address that is stated below for proper disposal

PHYWE Systeme GmbH & Co. KG
Department Customer Service
Robert-Bosch-Breite 10
D-37079 Göttingen

Phone +49 (0) 551 604-0
Fax +49 (0) 551 604-107

11 NOTES ON BATTERY AND RECHARGEABLE BATTERY DISPOSAL

As we sell batteries and rechargeable batteries or devices containing batteries and rechargeable batteries, we are obliged under the Battery Act (BattG) to inform you of the following: Batteries and rechargeable batteries may not be disposed of with household waste, but you are legally obliged to return used batteries and rechargeable batteries. Used batteries may contain harmful substances that can damage the environment or your health if they are not stored or disposed of properly. However, batteries also contain important raw materials such as iron, zinc, manganese or nickel and are recycled. You can either send the batteries back to us after use or return them free of charge in the immediate vicinity (e.g. in shops or at municipal collection centres). Batteries or rechargeable batteries that contain harmful substances are labelled with the symbol of a crossed-out dustbin.